



ECN - Tekniskt ändringsmeddelande

ECN-S04838

Datum: 2025-03-31


Utfärdare: Fredrik Nyblom

Sida 1 av 2

UPPDATERING ELBROMS REAL 9000 PLUS

Information om ändringar i Mercado Medics produkter eller tjänster.

ECR(s): S04838

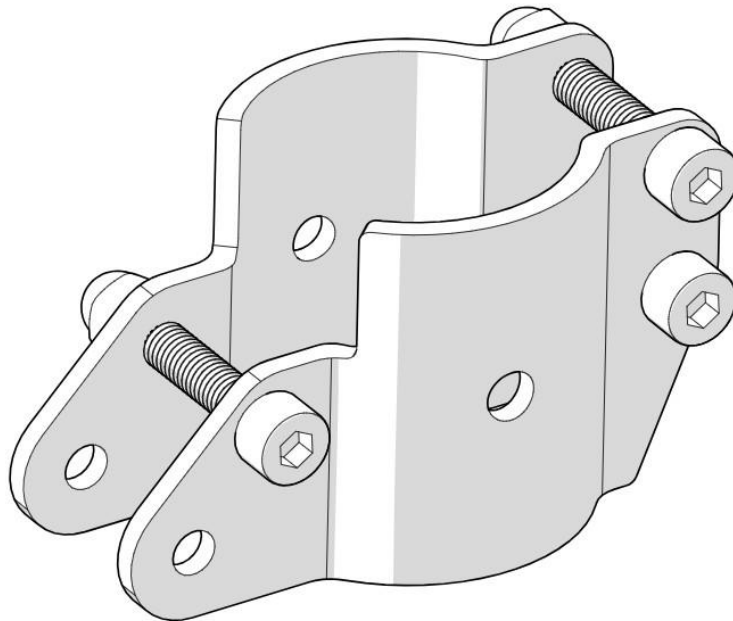

Sammanfattning

Elbroms Centrumfäste 9000f Eldon till REAL 9000 Plus har reviderats för enklare montering.

Artikelnummer som påverkas

802707	Elbroms Centrumfäste 9000f Eldon
804498	Eldonssats broms 24V
804482	Eldonssats broms Inkl Manbox utan tunga 24V
804548	Elbromskit kpl f manuell stol m manöverbox u tunga 24V
804621	Eldonssats broms Inkl manbox med tunga 24V
804627	Elbromskit kpl f manuell stol m manöverbox med tunga 24V
804870	Eldonssats broms f H-stativ Inkl Manbox utan tunga 24V

Bakgrund

Denna artikel har varit onödigt knepig att både tillverka och montera historiskt. Detta på grund av svetsmuttrarna som inte alltid satt rakt samt att formen på bockningen varit svår att få lika.

Detaljer kring ändringen

- Tillverkningsprocessen är ändrad till pressning i stället för kantpress.
- Svetsmuttrarna är ändrade till gängor och hattmuttrar.
- Tjockleken på plåten är ändrad från 3 till 4 mm.

Planerad tidpunkt för införande

2025-07-01

Kompatibilitet

Det går att byta till detta fäste på befintliga elbromsar.

Påverkan på kända risker

Ingen påverkan på befintliga risker. Inga nya risker uppstår.

Frågor

Vid frågor, vänligen kontakta vår produktchef Ronny Fogelqvist.

Kontaktuppgifter:

+46 708 27 96 14

ronny.fogelqvist@mercado.se



ECN - Engineering Change Notice

ECN-S04838

Date: 2025-03-31


Issued by: Fredrik Nyblom

Page 1 of 2

UPDATE OF THE ELECTRIC BRAKE REAL 9000 PLUS

Information about changes in Mercado Medic's products or services.

ECR(s): xxxx

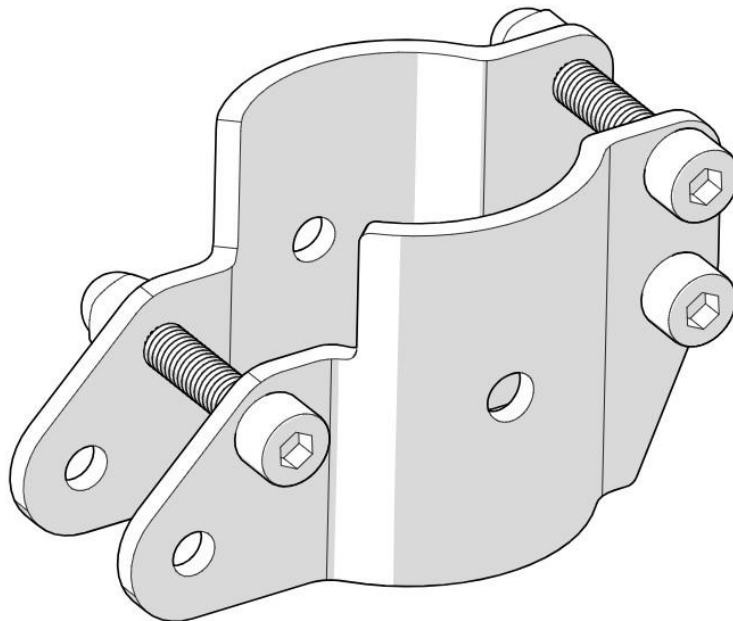

Summary

El. brake Centre attachment 9000 f electric kit for REAL 9000 Plus has been revised for easier assembly.

Affected article numbers

802707	El. brake Centre attachment 9000 f electric kit
804498	Actuator for el.brake 24V
804482	Actuator for el.brake incl.maneuverbox flat 24V
804548	El.brake kit for man.chair incl.maneuverbox flat 24V
804621	Actuator for el.brake incl. maneuverbox ext button 24V
804627	El.brake kit man.chair incl. maneuverbox ext button 24V
804870	Actuator f el brake f H-frame incl. maneuverbox flat 24V

Background

This item has historically been unnecessarily tricky to both produce and assemble. This is due to the weld nuts not always being properly aligned and the shape of the bending being difficult to make consistent.

Description of the change

- The manufacturing process has been changed to pressing instead of using a press brake.
- The weld nuts have been replaced with threads and cap nuts.
- The sheet metal thickness has been changed from 3 to 4 mm.

Estimated time for implementation

2025-07-01

Compatibility

It is possible to change to this bracket on existing electric brakes.

Safety information

No impact on existing risks. No new risks arise.

Questions

If you have any questions, please contact our Product Manager Ronny Fogelqvist.

Contact Information:

+46 708 27 96 14

ronny.fogelqvist@mercado.se